



Sterility Maintenance Products

Protect critical surfaces for aseptic processing

Protecting product contact surfaces from contamination during aseptic processing is critical to maintaining product sterility and ensuring patient safety. To address this need, STERIS offers a broad portfolio of sterility maintenance products, including sterilization wrapping, pre-sterilized cleanroom tools, specialty cleanroom apparel, and equipment covers.

Sterilization Wrapping Systems

The STERIS Sterilization Wrapping System is an optimal combination of autoclave bags, pouches, and covers designed specifically for your process. From the biggest, bulkiest stopper bowl to the smallest, most delicate filling needle, we will recommend the best wrapping configuration for your unique sterilization process. We offer proven designs that take into consideration parts wrapping, autoclave loading, sterile hold, transport, unwrapping, and aseptic presentation.

Pre-Sterilized Cleanroom Tools, Poly Bags, and Accessory Items

To simplify your cleanroom parts management and to reduce the need for autoclaving, STERIS provides a range of pre-sterilized items for your convenience. Whether you need a pair of sterile forceps to pick up vials, or a set of sterile mop bucket liners to reduce the frequency of autoclaving your mop buckets – STERIS can help design a solution that meets your needs. We also offer different packaging configurations to best facilitate the transfer of material through your facility and into the cleanroom.

Sterilization Bags

Sterilization Bags from STERIS are made from a flexible grade of Tyvek® nonwoven, spunbonded olefin material. They are available in a large range of sizes, and can be closed with an optional drawstring or with a heat sealer. Sterilization Bags are especially well suited for large, heavy, or irregularly shaped parts.



Sterilization Pouches with Film

Sterilization Pouches are comprised of Tyvek and film layers, forming a breathable pouch suitable for steam sterilization while offering a clear view of the contents inside. Options include a self-seal adhesive strip for closing the pouch and a chevron peel-open design for improved aseptic presentation.



Elasticized Bowl Covers | Bunker Covers | Hopper Covers

STERIS Elasticized Covers are a very simple yet highly effective way of protecting the product contact surfaces of stopper bowls, bunkers, and hoppers. The elasticized opening makes them easy to put on, easy to take off, and allows them to be left in place while the stopper bowl is installed on the fill line.



Elasticized Mini Covers

Easily cover small openings with STERIS Mini Covers. Made from a flexible Tyvek material with attached elastic, Mini Covers can quickly be manipulated onto or off of the open ends of process hoses, valves, fittings, tank ports, beakers, and carboys.



Equipment Covers

STERIS Equipment Covers can be used to cover tanks, racks, shelving, pallets, trolleys, carts, and other equipment in production and storage areas. Use of Equipment Covers provides an extra protective barrier between your equipment and the environment, and can help ensure GMP compliance.



Specialty Cleanroom Apparel

STERIS provides custom coveralls, bunny suits, and apparel packs/kits to meet the specific needs of your gowning procedures. Our line of specialty face masks include fog-prevention options and different fastener options.



Tyvek® is a registered trademark of DuPont or its affiliates.

Pre-Sterilized Cleanroom Tools

When tools or implements needed in the cleanroom cannot be autoclaved due to material compatibility or autoclave capacity issues, STERIS can provide you with pre-sterilized tools, such as pens, markers, forceps, scissors, and more. We can help you source a variety of items and provide the required packaging, sterilization, and documentation.





General Purpose Poly Bags

STERIS offers a range of poly bags to support the needs of cleanrooms for waste disposal, cleaning, transport, and storage. We provide a variety of trash bags, mop bucket liners, biohazard bags, zipper bags, and general purpose poly bags.



Sterility Maintenance Bags

Sterility Maintenance Bags provide a practical, cost-effective way to protect parts post-sterilization during transport and storage. STERIS Sterility Maintenance Bags are available in several different sizes, come with a self-seal adhesive strip, and can be easily opened without any cutting tools.

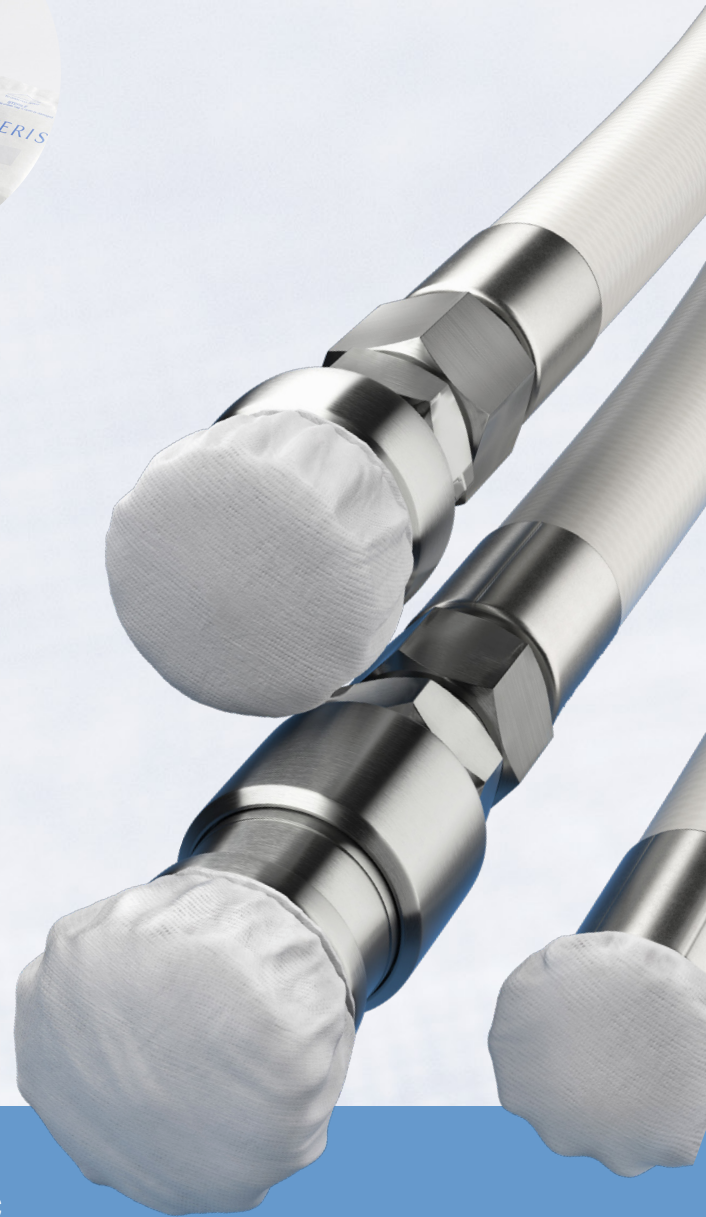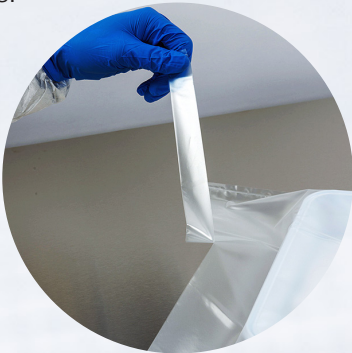

Selection of the optimal sterilization wrapping and other sterility maintenance products is a critical part of supporting your aseptic operations. With over 100 years of expertise in sterilization and contamination control, STERIS brings the industry's finest technical support and consultation to help you select and validate the right products for your contamination control program.